 SILVERGEAR

SMART SCALE PRO

BODY COMPOSITION & HEART RATE

User Manual

Thank you for purchasing our smart scale. This manual will provide instructions for effective and safe use of the product. Please read it carefully before use and store the manual in a safe place for future reference. If you ever pass the product on to others, this manual should be included.

1. Technical data

LCD display

Weight unit: kg/lb

Weight limit: 6-180 kg / 13.2-400 lb

Account number: no limitation (each account can add 24 users)

Baby mode: to weigh a baby

Size: 29.97 x 29.97 x 2.49 cm / 11.8 x 11.8 x 0.98 inch

Batteries: 4x 1.5V AA LR6 standard batteries

Support devices

- IOS 8.0 or higher
- Android 6.0 or higher

Recorded parameters

- Weight
- BMI
- Body fat
- Fat-free body weight
- Subcutaneous fat
- Heart rate
- Heart index
- Visceral fat
- Body water
- Skeletal muscle
- Muscle mass
- Bone mass
- Protein
- BMR
- Body age

Package content

- Smart scale
- Manual

2. Safety instructions

- The product is safe to use by children from the age of 10 and by persons with reduced physical or mental capabilities, or with lack of experience or knowledge, but only under supervision of an adult and with taking into account any hazards.
- Maintenance should be done by or under the supervision of an adult.
- Before use, inspect the device and its parts. Do not use in case of damage.
- Do not use to diagnose or treat any medical conditions.
- Do not use with medical implants such as pacemakers.
- Make sure the scale and your feet are dry before stepping on the scale, to avoid slipping.
- Do not stand on the edge of the scale or jump on it.
- Do not overload the scale (max 180 kg/400 lb).

- Place the product on a hard, flat surface.
- Do not drop the scale or drop objects on it, as this may damage the sensors.
- Do not immerse in water.
- Not recommended for pregnant women.
- Not recommended for infants, toddlers, and children under 10 years of age.
- Do not attempt to repair, disassemble or alter the product yourself.

Note:

- Any measurement obtained using this device, is for reference only and should not be considered as a medical opinion.
- Consult your GP or specialist before making changes to your diet, exercise plan or physical activities.
- To ensure the most accurate results, take your measurements at the same time of the day and under similar circumstances.

3. Preparations

Install Fitdays app

Search for 'Fitdays' in the Apple Store or in Google Play. Or scan this QR code: Download and install the app on your device. The app icon will appear on your phone or tablet after the installation is completed.



4. Operation

How to use Fitdays

Before using the smart scale, make sure you have entered all your personal data through the Fitdays app.

Turn on Bluetooth on your smart phone or tablet. Register your own ID by email or login with a social media account.

Add personal data and confirm OK.

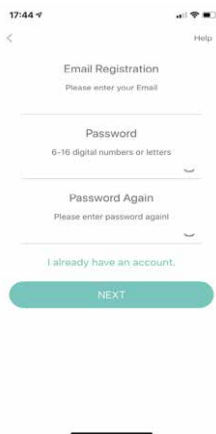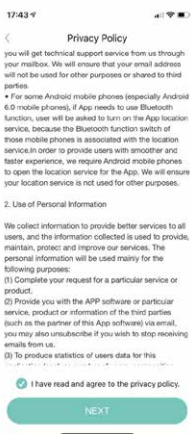Note: Entering the wrong gender, age, height, or corporeality level can result in inaccurate measurements.

17:45

17:49

 Male
  Female

Nickname

15075

Birthday

01/01/1995

Height

165cm

Target weight

56.1kg

Athlete mode

The personal information is only assessed as health data, Fitdays will keep you confidential

Confirm

 Male
  Female

Nickname

Fitdays

Birthday

29/11/1995

Height

163 cm


Target weight

54.7kg

Athlete mode

The personal information is only assessed as health data, Fitdays will keep you confidential

Confirm

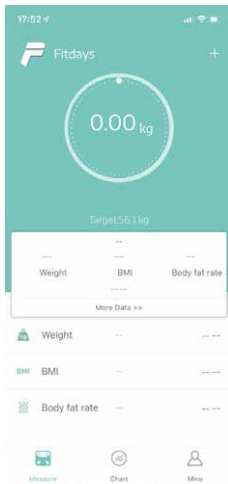Place the scale on a flat surface. Tap with your foot on the scale to activate the display. Wait until the display shows '0.0'(kg/lb). Step on the scale barefoot and with slightly damp feet for thorough measurements and full analysis. Make sure to touch all four metal squares. After the weight has been locked, all body composition data can be read in the Fitdays App by clicking More >> or Details.

Flat and hard surface

Clean and Bare feet

Carpet or uneven surface

With socks or shoes

Note: When stepping on the scale with socks on, only weight and BMI will be measured and displayed on the main page.

Usage tip:

Drink cold water when rising in the morning (ensure that the amount of water and its temperature are the same each day) and wait 30 minutes. Urinate and start the analysis with bare feet on the Silvergear scale.

Body balance

When the weight is locked, the scale will measure your body balance. Adjust your body to achieve equal weight distribution during tracking. The display flashes and displays your final weight, as well as your balance status.

28.0^L ≈ 24.6^R

Note: You should pay more attention to your balance, when the following data is shown:

- Weight L ≈ Weight R
- Health: Weight L > Weight R (or Weight L < Weight R)


Heart rate

After determining your body balance, the scale will measure the body's heart rate. The icon of a heart will flash when the scale is weighing and analysing. Keep standing on the scale until the measurement is finished and your phone has received all the data from the scale.

- You can find your personal data report on body weight, BMI, BMR, muscle, bone mass etc. by clicking on More >> or Details.
- When you click Charts, you can see your long-term tracking results.
- To delete results: click Charts to view your progress history and then click User Data. Use the calendar to select and slide to the left to delete the data.



Fitdays supports you in keeping a record for every function via the right icon. Fitdays can keep a track-record for recent, week, month and year. Share your results with friends and family Facebook, Instagram, or SMS.

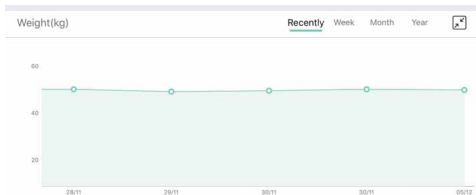Note: slide left or right or drag with two fingers to see more records.


- Add or switch users by clicking the icon upper-left, or click Me -> Manage users, then add personal data.
- The Smart Scale recognizes which user is on it and reminds him/her to switch to the correct user.
- Each account can support 24 users.

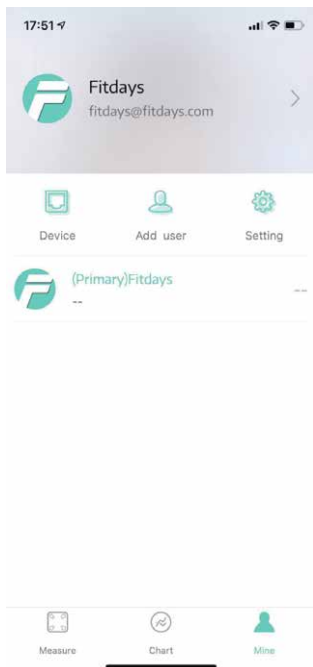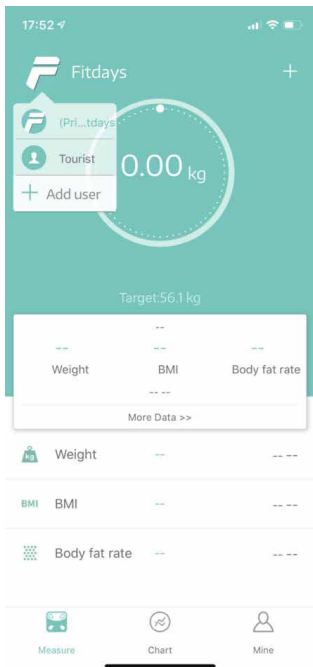Weighing a baby

Click + and Baby Mode. First, step on the scale yourself. When the scale has determined your weight, hold the baby to measure the baby's weight.

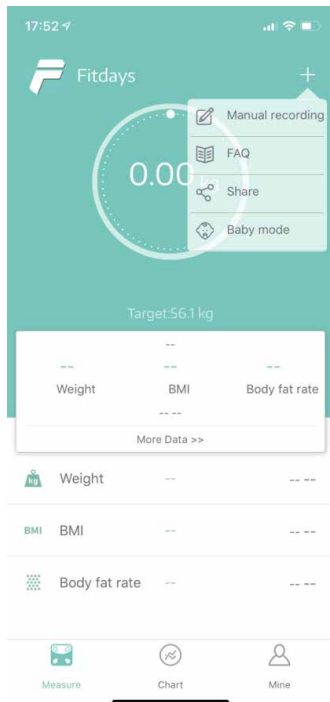17:59

Baby mode

58.50 kg

Weight of adult holding baby

Adult weighs first.

Please wait until the weight is steady before holding baby to weigh.

The APP automatically calculates the baby's weight.

17:57

Baby mode

6.80 kg

Baby weight

Please wait until the weight is steady before holding baby to weigh.

The APP automatically calculates the baby's weight.

5. Troubleshooting

1. Bluetooth is not connected

- Both the scale and smart phone must have **bluetooth** turned on.
- Check the **IOS** version of your smart phone; it must be **IOS 8.0** or higher / **Android 6.0** or higher.
- Click Me -> device find the scale

2. The scale is not measuring body fat data when weighing

- You have to be barefoot to measure body fat.
- The scale is not finished with weighing when the number displayed is flashing. Keep standing on the scale for about 10 seconds, until the number on the display stops flashing.
- **Bluetooth** is not turned on.
- **Bluetooth** is connected to a difference scale.

3. The scale is producing inconsistent measurements

- Check whether the scale is on a hard, flat surface.
- Check each sensor foot on the back to make sure nothing is stuck to the bottom.
- Remove the batteries to restart the scale.

4. The app is displaying an incorrect body fat rate

- The body parameter setting is incorrect, check to make sure that you enter the correct gender, height, and age.

5. How do I re-calibrate the scale (after moving it or changing the batteries)?

- Step on the scale to turn it on. Step off and let the number drop automatically to 0.0kg/lb to calibrate.

6. The scale does not turn on

- Remove the batteries and put them back. Make sure to insert the batteries as directed.

7. How to change units?

- You can change units via the Fitdays app. The scale will then show the same unit as the app.

1. Battery volume
2. Connected
3. Weight units
4. Body balance function
5. Body fat rate
6. Body water rate
7. Body Heart rate
8. Body weight
9. BMI
10. Muscle rate
11. Bone mass

No display

Functioning normally

Change battery

Lo

Change battery

EEr

Overload

[

Calibrating

6. Maintenance

- If necessary, clean the scale with a slightly damp towel and wipe it dry with a soft cloth.
- You can use cleaning alcohol to clean the sensors. Do not use it on other parts of the unit.
- Store the scale in a dry, cool place away from direct sunlight.

7. Disposal

Dispose of the article and the packaging materials in accordance with current local regulations.

This product is subject to the provisions of the European Directive 2012/19/EU.

The crossed out wheeled bin symbol means that this product shall not be disposed of with normal household waste. Electronic and Electrical Equipment not included in the selective sorting process are potentially dangerous for the environment and human health due to the presence of hazardous substances. Please dispose of responsibly at an approved waste or recycling facility.

The CE marking is a certification mark that indicates conformity with health, safety and environmental protection standards for products sold within the European Union. The CE marking is the manufacturer's declaration that the product meets the requirements of the applicable EC directives.

Recycling helps to reduce the consumption of resources and reduce the environmental impact. If you wish to discard the packaging, comply to the environmental legislation in your country.

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