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**Instruction manual**

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**Manual D'instruction**

 **ELEKTRISCHE PULSERENDE NEKMASSAGE**

 **ELECTRIC PULSE NECK MASSAGE**

 **PULS-NACKENMASSAGE**

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1. Don't use the product beside heart or on cephalosome, oral cavity, pudendum, skin disease position and other positions.
2. Patients with perceptual disturbance are prohibited to use.
3. Don't use this product in high temperature place, inflammable place, place with electromagnetic radiation and humid place.
4. Don't use when taking a shower, sweating and sleeping.

Please don't disassembly, repair and transform this product without permission to avoid the occurrence of failure or electric shock accident.

 Attention

1. Don't move the massager during the use procedure; firstly shut down the massager for massage position changing to avoid strong stimulus.
2. Please don't use it to children or people without expression consciousness ability.
3. Please stop use and consult the doctor when feeling uncomfortable in the process of using.
4. When you don't use it for a long time, please take out the batteries. Otherwise, it may cause battery leakage so that may lead to massager failure.

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# ELECTRIC PULSE NECK MASSAGE

## PRODUCT INTRODUCTION

QLT Neck Massager adopts the stream-line annular design in accordance with physiological curvature of human's cervical vertebra, uses the synchronous multichannel healthcare with trinity of low frequency electrotherapy, magnet therapy and annular traction to generate an efficient and composite energy field and to directly affect the deep tissue of organism and bone marrow. During the whole working process, it can always keep proper physiological position and posture and provide a scientific platform for the correct precaution and health care of cervical spondylosis. Provide you with easy, efficient and enjoyable multifunctional systematic health care at home and give you more economical enjoyment.

## SCOPE OF APPLICATION

The product can be applied to promote local blood circulation, relax local muscle and relieve fatigue. It is applicable to the health care and precaution of cervicodynia, numbness in shoulder and arms, cephalalgia and dizziness caused by cervical spondylosis.

## PRODUCT CHARACTERISTICS


1. According to ergonomic design, it is artistic and scientific and accords with neck curve;
2. Powered by two #7 batteries, safe and reliable;
3. Use synchronous multichannel healthcare with trinity of low frequency electrical impulse, magnetic effect and annular traction to generate an efficient and composite energy field;
4. Latest 3D intelligent anastomatic technology was adopted to regulate electrode easily in accordance with neck curve and provide more comfortable use;
5. 2 automatic composite massage modes and 3 manual massage modes can be selected according to diverse requirements to expand the scope of target users;
6. Multiple built-in health magnets can implement magnet therapy based on human acupoints, regulate the balance of yin-yang and promote microcirculation;
7. Configured with LED indicator light, its operation is simple and clear, and more humanized.

## PARTS INFORMATION

Mode Display Window

Conductive Pole Piece

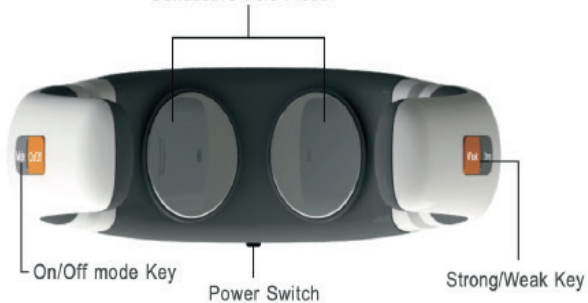On/Off mode Key

Power Switch

Strong/Weak Key

Battery  
Compartment

two #7 batteries

## OPERATION METHOD

### I. Power supply operation

Open battery cover and place the two #7 batteries into the battery compartment of massager based on the polarity.

### II. Functional operation

1) Press the power switch of massager to "ON", wipe the cervical vertebra and conductive pole piece with wet towel, and place the massager closely to neck. Please clean the sweat stain in cervical vertebra before use.

2) Press the "On/Off" key on the right hand shank of massager. The working procedure is started after hearing "beep" of the buzzer. Then indicator light 1 lights up. Press "Strong" key on the left hand shank once again and access into the setting automatic composite massage mode 1;

A. Automatic composite massage mode 1;

B. First gear is massage strength;

C. Timing: 15 minutes;

Press "On/off" key on the right hand shank to shut down and stop massage. Use of "MODE" key  
Mode switch key: six modes can be selected by pressing "MODE" on the right hand shank; the default boot setting is mode 1 (automatic composite mode 1) and the mode indicator light 1 lights up. Strength will be set as the first gear after every mode switch. The six massage modes are shown as below: Mode 1: automatic composite 1; mode indicator light 1 lights up;  
Mode 2: automatic composite 2; mode indicator light 2 lights up; Mode 3: manipulation; mode indicator light 3 lights up;

Mode 4: shiatsu; mode indicator light 1 and 3 light up at the same time;

Mode 5: malaxation; mode indicator light 2 and 3 light up at the same time; Mode 6: beating; mode indicator light 1, 2 and 3 light up at the same time;

Automatic composite mode 1 and 2 are manual modes behind the automatic cycle and will be switched automatically every minute. Automatic composite mode is recommended.

### Use of "STRONG" key

Strength increasing key: strength is set as 0-gear at each boot. Press "STRONG" key on the left hand shank once to increase one-gear strength. There are totally 16 gears. Please select the acceptable strength for massage. Low gear is recommended for first use.

### Use of "WEAK" key

Strength weakening key: press "WEAK" key on the left hand shank once to reduce one-gear strength. There are totally 16 gears. 10-20 minutes is recommended for daily use.

## NOTICE FOR LIFETIME AND REPLACEMENT OF BATTERY

1. AAA battery is available for the massager and can be used for about three months. (By one time for every day, 15 minutes for once).
2. Three continuous beeps mean shortage of power. Please replace a new battery.
3. When you don't see it for a long time, please take out the batteries. Otherwise, it may cause battery leakage so that may lead to massager failure.
4. Please don't use battery out of specified types.
5. Please pay attention to the positive and negative poles, distinguish the signs of "+" and "-".
6. Please shut off the power first to replace a new battery.
7. Please handle the used batteries according to relevant regulations of environment protection.

## CLEANING, MAINTENANCE AND STORAGE

1. AAA battery is available for the massager and can be used for about three months. (By one time for every day, 15 minutes for once).
2. Three continuous beeps mean shortage of power. Please replace a new battery.
3. When you don't see it for a long time, please take out the batteries. Otherwise, it may cause battery leakage so that may lead to massager failure.
4. Please don't use battery out of specified types.
5. Please pay attention to the positive and negative poles, distinguish the signs of "+" and "-".
6. Please shut off the power first to replace a new battery.
7. Please handle the used batteries according to relevant regulations of environment protection.

## SAFETY PRECAUTIONS

- The warning signs and figures in manual book are used to remind you to use this product safely and correctly and avoid the damage to yourself and others.
- Warning signs, figures and their meanings are shown as below:

Don't use it with other medical electronic instruments such as cardiac pacemaker, manual cardiopulmonary and other medical electronic instruments for life sustaining, electrocardiograph and others to avoid the occurrence of danger.


Using this product and high frequency surgical equipment at the same time may cause burn in the position of massager pole or damage the massager; using this product beside (1 meter) short wave or microwave therapeutic equipment may cause instable output of this product.

Pregnant women, women in menstrual period, skin allergy patients, patients with cardiopathy, dysrhythmia, malignancy, cerebrovascular disease, acute disease or other patients who are receiving medical treatment shall consult the doctor before using this product.

Warning signs	Meanings
 Dangerous	Incorrect use may cause the danger of casualties or serious injury.
 Warning	Incorrect use may cause the possibility of casualties or serious injury.
 Attention	Incorrect use may cause the possibility of casualties or goods damage.

### Figure signs

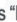Sign  means "PROHIBITION" (not allowed to do). The detailed prohibitive contents are presented by words or figures inside or beside . The meaning of left sign is "NO DISASSEMBLING".

Sign  means "ENFORCEMENT" (things must be abided by). The detailed enforced contents are presented by words or figures inside or beside . The meaning of left signal is "GENERAL ENFORCEMENT".

type BF Applied part.

The meaning of sign  is: attention, please refer to random file.

### Dangerous

Don't use it with other medical electronic instruments such as cardiac pacemaker, manual cardiopulmonary and other medical electronic instruments for life sustaining electrocardiograph and others to avoid the occurrence of danger.

Using this product and high frequency surgical equipment at the same time may cause burn in the position of massager pole or damage the massager; using this product beside (1 meter) short wave or microwave therapeutic equipment may cause instable output of this product.

### Warning

Pregnant women, women in menstrual period, skin allergy patients, patients with cardiopathy, dysarteriotony, malignancy, cerebrovascular disease, acute disease or other patients who are receiving medical treatment shall consult the doctor before using this product.

<ol style="list-style-type: none"> <li>1. Don't use the product beside heart or on cephalosome, oral cavity, pudendum, skin disease position and other positions.</li> <li>2. Patients with perceptual disturbance are prohibited to use.</li> <li>3. Don't use this product in high temperature place, inflammable place, place with electromagnetic radiation and humid place.</li> <li>4. Don't use when taking a shower, sweating and sleeping.</li> </ol>	
<p>Please don't disassembly, repair and transform this product without permission to avoid the occurrence of failure or electric shock accident.</p>	
 Attention	
<ol style="list-style-type: none"> <li>1. Don't move the massager during the use procedure; firstly shut down the massager for massage position changing to avoid strong stimulus.</li> <li>2. Please don't use it to children or people without expression consciousness ability.</li> <li>3. Please stop use and consult the doctor when feeling uncomfortable in the process of using.</li> <li>4. When you don't use it for a long time, please take out the batteries. Otherwise, it may cause battery leakage so that may lead to massager failure.</li> </ol>	

**Don't use it with other medical electronic instruments such as cardiac pacemaker, manual cardiopulmonary and other medical electronic instruments for life sustaining, electrocardiograph and others to avoid the occurrence of danger.**

**Using this product and high frequency surgical equipment at the same time may cause burn in the position of massager pole or damage the massager; using this product beside (1 meter) short wave or microwave therapeutic equipment may cause instable output of this product.**

**Pregnant women, women in menstrual period, skin allergy patients, patients with cardiopathy, dysarteriotony, malignancy, cerebrovascular disease, acute disease or other patients who are receiving medical treatment shall consult the doctor before using this product.**

### **IMPORTANT SAFETY INSTRUCTIONS**

the following basic safety guidance must always be followed when using electrical appliance to reduce the risk of fire, electric shock and/or injury to any person.

### **WARNING**

This device may be used by children aged eight and over, by people with reduced physical sensory or mental abilities, and by people lacking in experience and knowledge if they are appropriately supervised or given suitable instructions on the safe use of the device and understand the dangers its use entails.

Children should not be allowed to play with the device. Children should not clean or do maintenance on the device without supervision.

Consult your physician before using this product if:

- You are pregnant
- You have a pacemaker, an implant or other aids
- You suffer from diabetes
- You have a circulatory disorder
- You have concerns about your health

The elderly and infirm, should consult their doctor before using the product, even if there are no specific pre-existing conditions. Because skin friction can be painful, it is recommended that the device not be directly applied to the skin. Stop using immediately if you feel discomfort or pain.

Do not use this product at bedtime. The massager has a stimulating effect and can delay sleep. To avoid a short circuit, connect the device correctly and keep it away from metal (needles), waste and water. Use the device on a flat surface, if the device is not stable, it could tip over and cause injury.

To avoid injury and irreparable damage always sit while using the massager.

To avoid dangerous injury, never stand on the unit or use it forcefully.

Do not attempt to disassemble or repair the product. Repairs may only be carried out by suitably qualified technicians

- read all instructions fully and keep all instructions for future reference
- Not suitable for children
- Do not use while driving
- do not use if you are pregnant
- do not use if you have any of the following health conditions: heat sensitivity, circulation disorder, varicose vein, bruises, broken skin, skin irritation, muscle or tissue inflammation, vein inflammation, injury or any open wounds
- if you experience any pain of discomfort during the massage please stop immediately and consult your doctor
- Do not use if you have any of the following: Pacemaker, artificial joints, electronic implants
- Do not use the product if you have applied any oils or lotions to your skin
- Power off and disconnect the product from the mains after use
- Take precautions to make sure you do not fall asleep during use
- Don't carry this product by the power cord
- Never clean this product unless it is unplugged. Use a cloth dampened with a very small amount of water to clean. Can be used to wipe the product but it must be dried immediately. Do not plug this product back in until it is totally dry
- Do not operate under a blanket or pillow. Excessive heating can occur or cause fire, electric shock or injury to persons
- Never use this product near water or high humidity, never use this product unless you are totally dry
- Never immerse this product in or pour over with any liquid for any reason
- This appliance is suitable for personal household use only. Commercial or professional use will void the warranty. The supplier cannot be held responsible for injury or damage if the appliance has been used for anything other than its intended use
- This appliance must be kept out of the reach of children
- Always disconnect this appliance from the power supply when unattended.
- Do not leave this appliance unattended while plugged in
- Do not allow any leads of cords to rest in any position that might cause any trip or fall hazards
- Do not allow the power cord near any heated surfaces
- Do not use the appliance if it has been dropped or damaged in any way
- Check the power cord and plug regularly for any damage. If there is any damage to the power lead or plug it must be repaired by the manufacturer or a qualified professional
- Use only original parts if this appliance needs repair or maintenance
- If the appliance begins to malfunction in any way you must unplug and stop using it immediately
- Never try to inspect or open the product in any way
- Never operate or handle the appliance with wet hands
- Avoid using before bed - the massage can have a stimulating effect and may delay your sleep
- Keep all packaging provided with this appliance out of reach if children

## TECHNICAL DATA

Product specifications Name:

Neck Massager Model: 23C14/15000

Power supply: 3VDC 30mA (two #7 batteries)

Packing box size: 210×188×77mm Weight: about 410g

Normal operating conditions Environment temperature: +5°C~+40°C;

Relative humidity: ≤80%; Barometric pressure: 86kPa~106kPa;

Electrical performance

Pulse frequency: 1Hz~1000Hz.

Pulse width: 20us~400μ s. Pulse waveform: square wave.

Rated load impedance: 500Ω, allowing for ±10%;

Superficial magnetic field intensity: no less than 0.01T;

Timing period of massager: 15 minutes, tolerance allowed for ±10%;

4. Transportation and storage conditions Environment temperature: -20°C~+55°C

Relative humidity: ≤93%; Barometric pressure: 50kPa~106kPa.

Product features or design might change without notice

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