

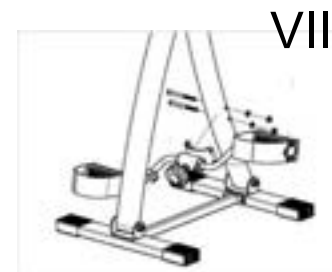
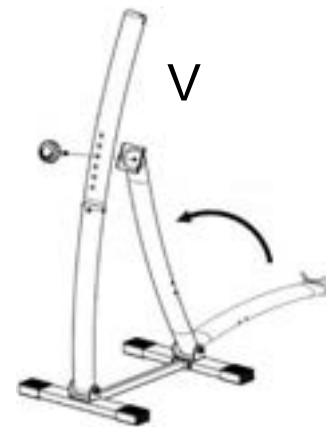
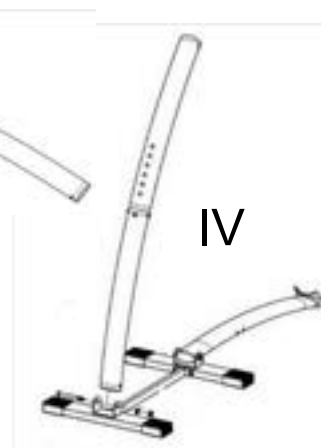
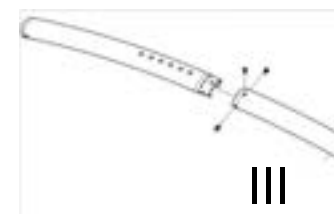
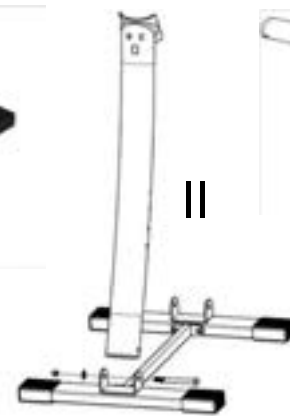
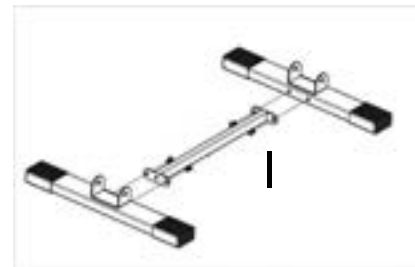
# BI-PEDALER®

-  *Instruction Manual*
-  *Manual de instrucciones*
-  *Mode d'emploi*
-  *Bedienungsanleitung*
-  *Istruzioni per l'uso*
-  *Manual de instruções*
-  *Gebruiksaanwijzing*

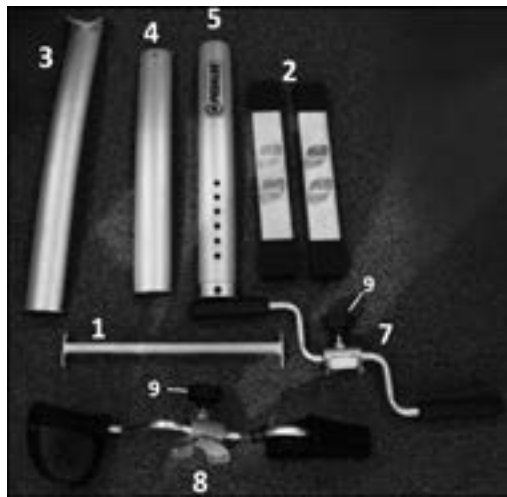International patents pending. All imitations will be prosecuted.  
BI-PEDALER® is a registered EU/CTM trade mark.
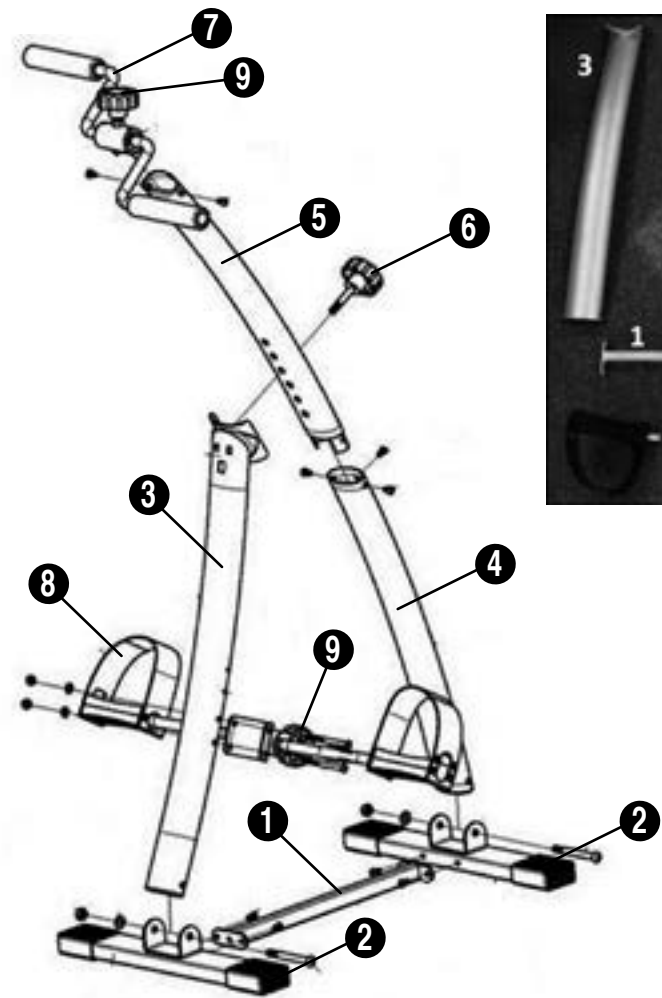
## BI-PEDALER® INSTRUCTION MANUAL PRODUCT COMPONENTS

### Product Description:

- |                               |                |
|-------------------------------|----------------|
| 1. Horizontal bar             | A. Screws x 3  |
| 2. Base tubes                 | B. Screws x 2  |
| 3. Rear bar                   | C. Washers x 2 |
| 4. Lower shaft                | D. Nuts x 2    |
| 5. Upper shaft                | E. Tool x 1    |
| 6. Height-adjustment knob     | F. Wrench x 1  |
| 7. Hand exerciser             | G. Screws x 2  |
| 8. Foot exerciser             | H. Washers x 2 |
| 9. Resistance-adjustment knob | I. Nuts x 2    |
|                               | J. Screws x 3  |

### Safety Warnings

Read the instruction manual before using BI-PEDALER® for the first time. If you have any doubts about your physical condition or have suffered from an illness, always consult your doctor before using BI-PEDALER®.

If you have a spinal disc herniation, consult your doctor before using the product.

If you have asthma or suffer from a respiratory disease, consult your doctor before using BI-PEDALER®. Always place the device on leveled base.

Ensure that children do not use BI-PEDALER® without adult supervision. Do not use if you are pregnant.

Do not drink alcoholic beverages while using the device. Drink water after use.

The maximum weight allowed on the BI-PEDALER® is 100 kg. For a safe, comfortable workout, leave a minimum of 0.6 m of free space

around the device.

Do not do any other exercises, since you may injure yourself by using the device incorrectly.

Warning: If any of the adjustment devices are left projecting, they could interfere with the user's movement.

### Class HC

### Set-up

First unpack all the parts.

### STEPS TO FOLLOW:

I. Attach the horizontal bar (1) to the 2 base tubes (2) using the A screws (which are stuck to the base tubes (2)). You will need tool E to attach them. Make sure they are firmly attached.

II. Attach the rear bar (3) to the base tube (2) as shown in the photo. Attach it with the B screws, C washers and D nuts using tool E and wrench F.

III. Attach the lower shaft (4) and the upper shaft (5) using the J screws, which are stuck to the bar. Use the E tool to attach them.

VI. Now attach the unit made up of the joined lower and upper shafts to the base tube (2), as shown in the photo, using the B screws, C washers and D nuts. You will need to use tool E and wrench F. Tighten enough to allow the bar to move on its axis.

V. Attach the joined lower shaft and upper shaft so it is resting on the rear bar (3), as shown in the photo. Attach them with the height adjustment knob (6). You can choose the height that works best for you from 6 different heights.

VI. Now take the hand exerciser (7) and remove the A screws. Place the hand exerciser in the hole in the shaft by inserting the tabs and attach it with the A screws using tool E. Make sure the exerciser is firmly attached.

VII. Place the foot exerciser (8) on the rear bar, as shown in the photo. Attach it using the 2 G screws, H washers and I nuts.


### Instructions for Use

BI-PEDALER® is not a professional piece of sports equipment. It is not recommended for use in a standing position. The user should be seated on a sofa or chair when using the device.

BI-PEDALER® has 6 different heights. Choose the one that works the best for you and insert the height-adjustment knob (6) in the corresponding hole.

The resistance of the hand and foot exercisers can be adjusted using the resistance-adjustment knob (9). Set the resistance that works best for you.

Important notice: After 7 minutes of use, do not touch the hand or foot exerciser to prevent burns.
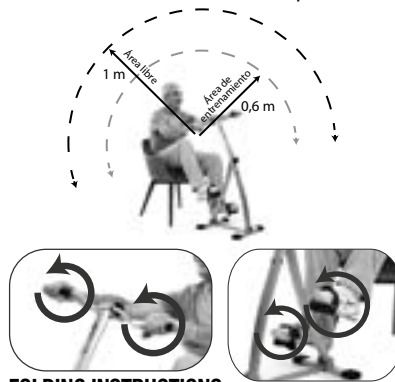
### FOLDING INSTRUCTIONS

1. Unscrew the resistance-adjustment knob (9) on the foot exerciser.
2. Unscrew the height-adjustment knob (6).
3. Fold the device as shown in Figure A.
4. Screw the height-adjustment knob (6) back into one of the holes.
5. Screw the resistance-adjustment knob (9) into the top hole on the shaft. See Figure B.

To assemble BI-PEDALER® again, repeat the steps in reverse order.

BI-PEDALER® comes with a magnetic bracelet as a free gift. For maintenance purposes, keep it from becoming wet.

### Precautions:

- This product is not a toy and should be kept out of reach of children.
- If you are pregnant or think you might be, do not use the device or use it only after consulting your doctor.
- Patients with heart problems or implants should not use the device.
- Patients with haemophilia, thrombocytopenia or other blood disorders should not use the device unless authorized by their doctor.
- The permanent magnet may cause allergies. If this occurs, stop wearing the bracelet.

### Quality Warranty:

This product is covered by a warranty against manufacturing defects subject to the time periods stipulated by the legislation in force in each country.

This warranty does not cover damages resulting from inadequate use, negligent commercial use, abnormal wear and tear, accidents or improper handling.

### Australia & New Zealand only

This product has the benefit of certain Consumer Guarantees. These are prescribed by the Australian Consumer Law & the New Zealand Consumer Guarantees Act 1993 both of which provide protection for consumers. There is no express warranty for this product in Australia or New Zealand. The above paragraph refers to other countries.

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ISL GmbH
Best Direct [International] Benelux B.V.
Best Direct International Ltd,
ISL ITALY S.r.l.